


Ascent Unboxing and Setup




The Ascent instrument is heavy. Be careful when lifting.


1. Remove the accessory tray, guides, and instrument from the box and place them on a stable, level bench out of direct sunlight.

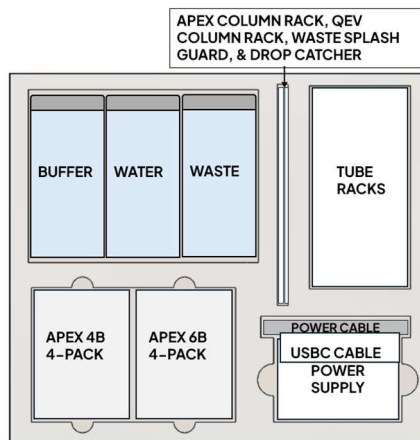

2. Remove the red shipping tape holding the carriage and tube rack platform in place.

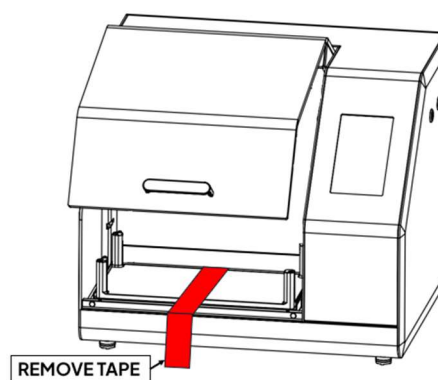

3. Connect each bottle to the corresponding tube on the side of the instrument. Firmly press the tubing into the fitting until it stops. Place the waste splash guard on the waste.

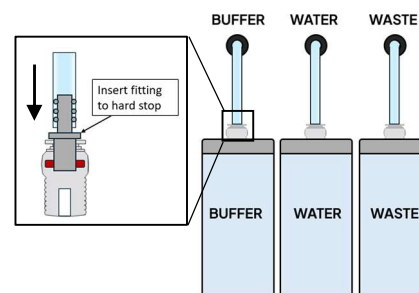

-
4. Plug the power supply into the back of the instrument and turn on the power switch next to the power input. Refer to the Ascent quick user guide to begin using the instrument.
-

